




PATIENT

Routine X-Ray Kit

INF-17, INF-17-SCB, INF17-UF



INF-17 Open Cell Foam

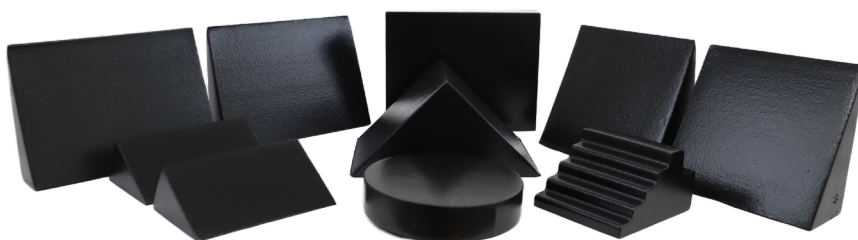

INF-17-SCB Open Cell Foam with ScanCoat®

Each 10-Piece Kit Includes:

- (1) 112 - Circular Disc, 1.5”H x 7” Diameter
- (1) 114 - Rectangle, 2”H x 7.5”W x 9.5”L
- (2) 116 - 34 degree Wedge, 3”H x 4”W x 7”L
- (2) 117 - 16 degree Wedge, 2.5”H x 7”W x 10”L
- (1) 118 - 45 degree Wedge, 7”H x 7”W x 7”L
- (1) 121 - Oblique Finger Block, 3.13”H x 5.5”W x 5”L
- (2) 122 - 22 degree Wedge, 3.5”H x 7”W x 8”L

Routine X-Ray Kit Available Finishes

- Open Cell Foam (Part #INF-17)
Uncovered standard foam, charcoal colored, open cell foam, radiolucent, fire retardant, economic and useful for general positioning.
- Flexible ScanCoat (Part #INF-17-SCB)
State of the art sprayed vinyl coating, withstands tearing or peeling, skid resistant, resists fluid absorption and very low attenuation. This vinyl coating is ideal for imaging compatibility, cleanability, and contamination control.
- Closed Cell Foam (Part #INF-17-UF)
Ultrafoam: Firm, closed cell foam, charcoal colored, resists fluid absorption/contamination and is fire retardant. Ultrafoam construction is ideal for firm support during positioning applications (17-UF).

*Please note: Covered items (-8, -20 or -SCB) are not recommended for low KVP techniques as the edges/seams may show. Items composed of Ultrafoam (-UF) that have a height higher than 3 inches will have a glue seam that may attenuate.

Elevate your imaging capabilities and patient care standards with INFAB’s premium X-Ray Patient Positioning 10-piece Kit for enhanced precision and efficiency. Our routine x-ray positioning kit includes an assortment of foam positioning aids commonly used in a variety of imaging applications including general positioning, CT scans, MRI, nuclear medicine or ultrasound.

- Versatile routine x-ray positioning kit
- Includes 10 lightweight foam positioning aids
- Ensures proper patient positioning and imaging efficiency
- Created specifically to meet the patient positioning needs of your imaging department
- Patient Positioning Devices are available uncoated and coated
- All coverings/material are Latex Free, MRI Safe & Biocompatible

| | |
|--------------------|---|
| Model | INF-17 Open Cell Foam INF-17-SCB Open Cell Foam with ScanCoat® INF-17-UF Closed Cell Foam |
| Description | Routine X-ray Kit |
| Modalities | CT Scan, General Positioning, MRI, Nuclear Medicine, Ultrasound |
| Sterile | N |
| Latex Free | Y |
| Single Use | N |



REV002

1040 Avenida Acaso
Camarillo, CA 93012, USA

tel: 805.987.5255
fax: 805.482.8424

info@infabcorp.com
www.infabcorp.com

